

# Bellman® Audio

## NeckLoop BE9128



## For listening with hearing aids

Light and discrete neck loop with high signal output

For high quality and easy listening with your hearing aids. Works with any sound source that features a 3.5 mm headphones output. Simply plug the Bellman® Audio Neck loop into your MP3 player, mobile computer, hands free or personal amplifier and put your hearing aid in "T" position.

The Bellman® Audio Neck Loop is extremely discrete and light weight, and can easily be worn under your shirt.

### Facts Bellman® Audio Neck Loop

- Neck strap length: 70 cm
- Cable length: 90 cm
- Connector: 3.5 mm tele plug (stereo), 90 degrees angle
- Colour: Black
- Weight: 30 g
- Size: width 3,5 cm, height 4,0 cm, depth 1,0 cm
- Features a safety breakaway neck strap
- Impedance: 2 x 5 ohm
- Output: 100 mA/m at 10 cm distance and 2 x 1.0 Vrms incoming signal

**Bellman® audio**  
FOR EASY COMMUNICATION



### Features

**Ultra light**  
and discrete

**Dual channel**  
for high signal output

### Sound sources



MP3 Player



Mobile Phone




Mobile Computer



Personal Amplifier

### How to connect

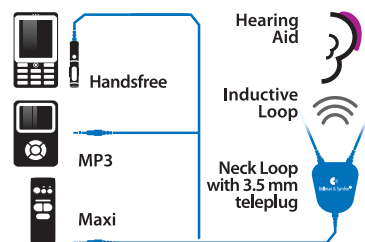

**Bellman & Symfon®**

Södra Långebergsgatan 30, 421 32 Västra Frölunda  
Tel +46 (0)31 68 28 20, Fax +46 (0)31 68 28 90  
www.bellman.se, info@bellman.se